


---

**Instructions for use**  
**2-CAT (N-D) Research ELISA™**

**REF**

**BA E-5500**

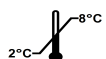

**RUO**

For research  
use only –  
Not for use  
in diagnostic  
procedures

## **2-CAT (N-D) Research ELISA**

### **1. Intended use and principle of the test**

Enzyme Immunoassay for the quantitative determination of noradrenaline (norepinephrine) and dopamine. Flexible test system for various biological sample types and volumes.

Noradrenaline (norepinephrine) and dopamine are extracted by using a cis-diol-specific affinity gel, acylated and then converted enzymatically.

The competitive ELISA kit uses the microtiter plate format. The antigen is bound to the solid phase of the microtiter plate. The derivatized standards, controls and samples and the solid phase bound analytes compete for a fixed number of antibody binding sites. After the system is in equilibrium, free antigen and free antigen-antibody complexes are removed by washing. The antibody bound to the solid phase is detected by an anti-rabbit IgG-peroxidase conjugate using TMB as a substrate. The reaction is monitored at 450 nm.

Quantification of unknown samples is achieved by comparing their absorbance with a standard curve prepared with known standard concentrations.

### **2. Procedural Cautions, Guidelines and Warnings**


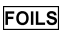
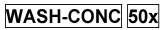
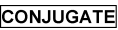

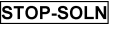








- (1) This kit is intended for professional use only. Users should have a thorough understanding of this protocol for the successful use of this kit. Only the test instruction provided with the kit is valid and has to be used to run the assay. Reliable performance will only be attained by strict and careful adherence to the instructions provided.
- (2) The principles of Good Laboratory Practice (GLP) have to be followed.
- (3) In order to reduce exposure to potentially harmful substances, wear lab coats, disposable latex gloves and protective glasses where necessary.
- (4) All kit reagents and specimens should be brought to room temperature and mixed gently but thoroughly before use. Avoid repeated freezing and thawing of reagents and specimens.
- (5) For dilution or reconstitution purposes, use deionized, distilled, or ultra-pure water.
- (6) The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch with desiccant and used in the frame provided. Microtiter strips which are removed from the frame for usage should be marked accordingly to avoid any mix-up.
- (7) Duplicate determination of sample is highly recommended to be able to identify potential pipetting errors.
- (8) Once the test has been started, all steps should be completed without interruption. Make sure that the required reagents, materials and devices are prepared ready at the appropriate time.
- (9) Incubation times do influence the results. All wells should be handled in the same order and time intervals.
- (10) To avoid cross-contamination of reagents, use new disposable pipette tips for dispensing each reagent, sample, standard and control.
- (11) A standard curve must be established for each run.
- (12) The controls should be included in each run and fall within established confidence limits. The confidence limits are listed in the QC-Report.
- (13) Do not mix kit components with different lot numbers within a test and do not use reagents beyond expiry date as shown on the kit labels.
- (14) Avoid contact with Stop Solution containing 0.25 M H<sub>2</sub>SO<sub>4</sub>. It may cause skin irritation and burns. In case of contact with eyes or skin, rinse off immediately with water.
- (15) TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them.
- (16) For information on hazardous substances included in the kit please refer to Safety Data Sheet (SDS). The Safety Data Sheet for this product is made available directly on the website of the manufacturer or upon request.
- (17) Kit reagents must be regarded as hazardous waste and disposed according to national regulations.

### **3. Storage and stability**

Store the unopened reagents at 2 - 8 °C until expiration date. Do not use components beyond the expiry date indicated on the kit labels. Once opened the reagents are stable for 1 month when stored at 2 - 8 °C. Once the resealable pouch has been opened, care should be taken to close it tightly with desiccant again.

## 4. Materials

### 4.1 Content of the kit

<b>BA D-0032</b>	 96	<b>Microtiter Plate</b> - Ready to use
Content:	1 x 96 wells, empty in a resealable pouch	
<b>BA D-0090</b>	 FOILS	<b>Adhesive Foil</b> - Ready to use
Content:	Adhesive Foils in a resealable pouch	
Volume:	2 x 4 foils	
<b>BA E-0030</b>	 WASH-CONC 50x	<b>Wash Buffer Concentrate</b> - Concentrated 50x
Content:	Buffer with a non-ionic detergent and physiological pH	
Volume:	2 x 20 ml/vial, light purple cap	
<b>BA E-0040</b>	 CONJUGATE	<b>Enzyme Conjugate</b> - Ready to use
Content:	Goat anti-rabbit immunoglobulins, conjugated with peroxidase	
Volume:	2 x 12 ml/vial, red cap	
<b>BA E-0055</b>	 SUBSTRATE	<b>Substrate</b> - Ready to use
Content:	Chromogenic substrate containing tetramethylbenzidine, substrate buffer and hydrogen peroxide	
Volume:	2 x 12 ml/black vial, black cap	
<b>BA E-0080</b>	 STOP-SOLN	<b>Stop Solution</b> - Ready to use
Content:	0.25 M sulfuric acid	
Volume:	2 x 12 ml/vial, light grey cap	
Hazards identification:		
	H290 May be corrosive to metals. H314 Causes severe skin burns and eye damage.	
<b>BA E-0231</b>	 NAD NMN	<b>Noradrenaline Microtiter Strips</b> - Ready to use
Content:	1 x 96 well (12x8) antigen precoated microwell plate in a resealable yellow pouch with desiccant	
<b>BA E-0331</b>	 DOP	<b>Dopamine Microtiter Strips</b> - Ready to use
Content:	1 x 96 well (12x8) antigen precoated microwell plate in a resealable green pouch with desiccant	
<b>BA E-5210</b>	 NAD-AS	<b>Noradrenaline Antiserum</b> - Ready to use
Content:	Rabbit anti-noradrenaline antibody, yellow coloured	
Volume:	1 x 6 ml/vial, yellow cap	
<b>BA E-5310</b>	 DOP-AS	<b>Dopamine Antiserum</b> - Ready to use
Content:	Rabbit anti-dopamine antibody, green coloured	
Volume:	1 x 6 ml/vial, dark green cap	
<b>BA R-0050</b>	 ADJUST-BUFF	<b>Adjustment Buffer</b> - Ready to use
Content:	TRIS buffer	
Volume:	1 x 4 ml/vial, green cap	
<b>BA R-4617</b>	 TE-BUFF	<b>TE Buffer</b> - Ready to use
Content:	TRIS-EDTA buffer	
Volume:	1 x 4 ml/vial, brown cap	
<b>BA R-6618</b>	 EXTRACT-PLATE 48	<b>Extraction Plate</b> - Ready to use
Content:	2 x 48 well plates coated with boronate affinity gel in a resealable pouch	

**BA R-6619** HCL **Hydrochloric Acid** - Ready to use  
 Content: 0.025 M Hydrochloric Acid, yellow coloured  
 Volume: 1 x 20 ml/vial, dark green cap

**Standards and Controls** - Ready to use

Cat. no.	Component	Colour/ Cap	Concentration ng/ml		Concentration nmol/l		Volume/ Vial
			NAD	DOP	NAD	DOP	
<b>BA R-5601</b>	<span style="border: 1px solid black; padding: 0 2px;">STANDARD A</span>	white	0	0	0	0	4 ml
<b>BA R-5602</b>	<span style="border: 1px solid black; padding: 0 2px;">STANDARD B</span>	light yellow	0.2	0.5	1.2	3.3	4 ml
<b>BA R-5603</b>	<span style="border: 1px solid black; padding: 0 2px;">STANDARD C</span>	orange	0.6	1.5	3.5	9.8	4 ml
<b>BA R-5604</b>	<span style="border: 1px solid black; padding: 0 2px;">STANDARD D</span>	dark blue	2	5	12	33	4 ml
<b>BA R-5605</b>	<span style="border: 1px solid black; padding: 0 2px;">STANDARD E</span>	light grey	8	20	47	131	4 ml
<b>BA R-5606</b>	<span style="border: 1px solid black; padding: 0 2px;">STANDARD F</span>	black	32	80	189	522	4 ml
<b>BA R-5651</b>	<span style="border: 1px solid black; padding: 0 2px;">CONTROL 1</span>	light green	Refer to QC-Report for expected value and acceptable range!				4 ml
<b>BA R-5652</b>	<span style="border: 1px solid black; padding: 0 2px;">CONTROL 2</span>	dark red					4 ml

Conversion: Noradrenaline (ng/ml) x 5.91 = Noradrenaline (nmol/l)  
 Dopamine (ng/ml) x 6.53 = Dopamine (nmol/l)

Content: Acidic buffer with non-mercury stabilizer, spiked with defined quantity of noradrenaline and dopamine

**BA R-6611** ACYL-BUFF **Acylation Buffer** - Ready to use  
 Content: Buffer with light alkaline pH for the acylation  
 Volume: 1 x 20 ml/vial, white cap

**BA R-6612** ACYL-REAG **Acylation Reagent** - Ready to use  
 Content: Acylation reagent in DMF and DMSO  
 Volume: 1 x 3 ml/vial, light red cap

Hazards identification:



H226 Flammable liquid and vapour.  
 H360D May damage the unborn child.  
 H312 + H332 Harmful in contact with skin or if inhaled.  
 H319 Causes serious eye irritation.

**BA R-6614** COENZYME **Coenzyme** - Ready to use  
 Content: S-adenosyl-L-methionine  
 Volume: 1 x 4 ml/vial, purple cap

**BA R-6615** ENZYME **Enzyme** - Lyophilized  
 Content: Catechol-O-methyltransferase  
 Volume: 4 vials, pink cap

#### 4.2 Additional materials and equipment required but not provided in the kit

- Calibrated precision pipettes to dispense volumes between 1 – 750 µl; 1 ml
- Microtiter plate washing device (manual, semi-automated or automated)
- ELISA reader capable of reading absorbance at 450 nm and if possible 620 – 650 nm
- Shaker (shaking amplitude 3 mm; approx. 600 rpm)
- Temperature controlled incubator (37 °C) or similar heating device
- Absorbent material (paper towel)
- Water (deionized, distilled, or ultra-pure)
- Vortex mixer

#### 5. Sample collection and storage

Storage: up to 6 hours at 2 – 8 °C; for longer periods (up to 6 months) at -20 °C or -80 °C.

*Advice for the preservation of the biological sample:* to prevent catecholamine degradation, add EDTA (final concentration 1 mM) and sodium metabisulfite (final concentration 4 mM) to the sample.

## 6. Test procedure

Allow all reagents and samples to reach room temperature and mix thoroughly by gentle inversion before use. Duplicate measurements are recommended. It is recommended to number the strips of the microwell plate before usage to avoid any mix-up.

The binding of the antiserum and the enzyme conjugate and the activity of the enzyme are temperature dependent, and the absorbance may vary if a thermostat is not used. The higher the temperature, the higher the absorbance will be. Varying incubation times will have a similar influence on the absorbance. The optimal temperature during the Enzyme Immunoassay is between 20 – 25 °C.

■ *In case of overflow, read the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to 405 nm*

### 6.1 Preparation of reagents

#### Wash Buffer

Dilute the 20 ml Wash Buffer Concentrate with water (deionized, distilled, or ultra-pure) to a final volume of 1000 ml.

Storage: 1 month at 2 – 8 °C

#### Enzyme Solution

Reconstitute the content of the vial labelled 'Enzyme' with 1 ml water (deionized, distilled, or ultra-pure) and mix thoroughly. Add 0.3 ml of Coenzyme followed by 0.7 ml of Adjustment Buffer. The total volume of the Enzyme Solution is 2.0 ml.

■ *The Enzyme Solution has to be prepared freshly prior to the assay (not longer than 10 - 15 minutes in advance). Discard after use!*

#### Noradrenaline Microtiter Strips and Dopamine Microtiter Strips

In rare cases residues of the blocking and stabilizing reagent can be seen in the wells as small, white dots or lines. These residues do not influence the quality of the product.

### 6.2 Sample preparation

The 2-CAT (N-D) Research ELISA is a flexible test system for various biological sample types and volumes. It is not possible to give a general advice how to prepare the samples. However, the following basics should help the researcher to fit the protocol to his specific needs.

- Avoid excess of acid: excess of acid might exceed the buffer capacity of the extraction buffer. A pH > 7.0 during the extraction is mandatory.
- Prevent catecholamine degradation by adding preservatives to the sample (please refer to Sample collection and storage).
- Avoid chaotropic chemicals like perchloric acid. The high salt content might reduce the recovery of catecholamines. If your samples already contain high amounts of perchloric acid, neutralize the sample prior to the extraction step.
- Tissue samples can be homogenised in 0.01 N HCl in the presence of EDTA and sodium metabisulfite. Under these conditions, catecholamines are positively charged which reduces binding to proteins and optimizes solubility.
- Avoid samples that contain substances with a cis-diol structure. These will reduce the recovery of the catecholamines.
- It is advisable to perform a "Proof of Principle" to determine the recovery of the catecholamines in your samples. Prepare a stock solution of noradrenaline and dopamine. Add small amounts (to change the native sample matrix as less as possible) of the stock solutions to the sample matrix and check the recovery.
- The used sample volume determines the sensitivity of the test. Determine the sample volume needed to determine the catecholamines in your sample by testing different amounts of sample volume.

*If you need any support in establishing a protocol for your specific purposes, do not hesitate to contact the manufacturer directly!*

### 6.3 Extraction and Acylation

The 2-CAT (N-D) Research ELISA offers a flexible test system for various biological sample types and volumes. Step 1 of the extraction procedure depends on the sample volume:

- in case you have sample volumes between 1 – 100 µl follow **1.1**
- in case you have sample volumes between 100 – 500 µl follow **1.2**
- in case you have sample volumes between 500 – 750 µl follow **1.3**

■ ***Within a run it is only possible to measure samples with the same volume***

1.	1.1 Sample volume 1 – 100 µl	1.2 Sample volume 100 – 500 µl	1.3 Sample volume 500 – 750 µl
	Pipette into the respective wells of the Extraction Plate: <b>20 µl standards, 20 µl controls and 1 – 100 µl sample.</b> Fill up each well with water (deionized, distilled, or ultra-pure) to a <b>final volume</b> of 100 µl [e.g. 20 µl standard plus 80 µl water (deionized, distilled, or ultra-pure)].	Pipette into the respective wells of the Extraction Plate: <b>20 µl standards, 20 µl controls and 100 – 500 µl sample.</b> Fill up each well with water (deionized, distilled, or ultra-pure) to a <b>final volume</b> of 500 µl [e.g. 20 µl standard plus 480 µl water (deionized, distilled, or ultra-pure)].	Pipette into the respective wells of the Extraction Plate: <b>20 µl standards, 20 µl controls and 500 – 750 µl sample.</b> Fill up each well with water (deionized, distilled, or ultra-pure) to a <b>final volume</b> of 750 µl [e.g. 20 µl standard plus 730 µl water (deionized, distilled, or ultra-pure)].
2.	Pipette <b>25 µl</b> of <b>TE Buffer</b> into all wells		
3.	Cover the plate with <b>Adhesive Foil</b> . Shake <b>60 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).		
4.	Remove the foil and empty the plate. Blot dry by tapping the inverted plate on absorbent material.		
5.	Pipette <b>1 ml</b> of <b>Wash Buffer</b> into all wells.		
6.	Shake <b>5 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).		
7.	Blot dry by tapping the inverted plate on absorbent material.		
8.	<b>Wash one more time</b> as described (step 5, 6 and 7)!		
9.	Pipette <b>150 µl</b> of <b>Acylation Buffer</b> into all wells.		
10.	Pipette <b>25 µl</b> of <b>Acylation Reagent</b> into all wells.		
11.	Shake <b>20 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).		
12.	Empty the plate and blot dry by tapping the inverted plate on absorbent material.		
13.	Pipette <b>1 ml</b> of <b>Wash Buffer</b> into all wells.		
14.	Shake <b>5 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).		
15.	Blot dry by tapping the inverted plate on absorbent material.		
16.	<b>Wash one more time</b> as described (step 13, 14, 15).		
17.	Pipette <b>150 µl</b> of <b>Hydrochloric Acid</b> into all wells.		
18.	Cover plate with <b>Adhesive Foil</b> . Shake <b>10 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).		
	<b>Do not decant the supernatant thereafter!</b>		
	<b>140 µl of the supernatant is needed for the subsequent enzymatic conversion</b>		

#### 6.4 Enzymatic Conversion

1.	Pipette <b>140 µl</b> of the <b>extracted standards, controls</b> and <b>samples</b> into the respective wells of the <b>Microtiter Plate</b> .				
2.	Add <b>50 µl</b> of <b>Enzyme Solution</b> (refer to 6.1) to all wells.				
3.	Cover plate with <b>Adhesive Foil</b> . Shake <b>1 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).				
4.	Incubate for <b>2 h</b> at <b>37°C</b> . The following volumes of the supernatants are needed for the subsequent ELISA:				
	<table border="1" style="display: inline-table; margin-right: 20px;"> <tr> <td style="padding: 2px 10px;"><b>Noradrenaline</b></td> <td style="padding: 2px 10px;"><b>90 µl</b></td> </tr> </table> <table border="1" style="display: inline-table;"> <tr> <td style="padding: 2px 10px;"><b>Dopamine</b></td> <td style="padding: 2px 10px;"><b>90 µl</b></td> </tr> </table>	<b>Noradrenaline</b>	<b>90 µl</b>	<b>Dopamine</b>	<b>90 µl</b>
<b>Noradrenaline</b>	<b>90 µl</b>				
<b>Dopamine</b>	<b>90 µl</b>				

## 6.5 Noradrenaline and Dopamine ELISA

1. Pipette <b>90 µl</b> of <b>standards, controls</b> and <b>samples</b> from the <b>Enzyme Plate</b> (refer to 6.4) into the respective pre-coated <b>Microtiter Strips</b> (* <sup>1</sup> ).
2. Pipette <b>50 µl</b> of the respective <b>Antiserum</b> (* <sup>2</sup> ) into all wells.
3. Cover the plate with <b>Adhesive Foil</b> . Shake <b>1 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).
4. Incubate for <b>15 – 20 h</b> (overnight) at <b>2 – 8 °C</b> .
5. Remove the foil. Discard or aspirate the content of the wells. Wash the plate <b>4 x</b> by adding <b>300 µl</b> of <b>Wash Buffer</b> , <b>discarding</b> the content and <b>blotting dry each time</b> by tapping the inverted plate on absorbent material.
6. Pipette <b>100 µl</b> of <b>Enzyme Conjugate</b> into all wells.
7. Cover the plate with <b>Adhesive Foil</b> . Incubate <b>30 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm).
8. Remove the foil. Discard or aspirate the content of the wells. Wash the plate <b>4 x</b> by adding <b>300 µl</b> of <b>Wash Buffer</b> , <b>discarding</b> the content and <b>blotting dry each time</b> by tapping the inverted plate on absorbent material.
9. Pipette <b>100 µl</b> of <b>Substrate</b> into all wells.
10. Incubate <b>20 - 30 min</b> at <b>RT</b> (20 – 25 °C) on a <b>shaker</b> (approx. 600 rpm). ■ <b>Avoid exposure to direct sunlight!</b>
11. Pipette <b>100 µl</b> of <b>Stop Solution</b> into all wells.
12. <b>Read</b> the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to <b>450 nm</b> (if available a reference wavelength between 620 nm and 650 nm is recommended).

▲ (\*<sup>1</sup>): **Noradrenaline Microtiter Strips, Dopamine Microtiter Strips**

▲ (\*<sup>2</sup>): **Noradrenaline Antiserum, Dopamine Antiserum**

## 7. Calculation of results

The standard curve from which the concentrations of the samples can be read off, is obtained by plotting the absorbance readings (calculate the mean absorbance) measured for the standards (linear, y-axis) against the corresponding standard concentrations (logarithmic, x-axis).

Use a non-linear regression for curve fitting (e.g. spline, 4- parameter, akima).

■ *This assay is a competitive assay. This means: the OD-values are decreasing with increasing concentrations of the analyte. OD-values found below the standard curve correspond to high concentrations of the analyte in the sample and have to be reported as being positive.*

■ **The concentrations of the samples taken from the standard curve have to be multiplied by a correction factor.**

$$\text{Correction factor} = \frac{20 \mu\text{l (volume of standards extracted)}}{\text{sample volume } (\mu\text{l}) \text{ extracted}}$$

### Example:

750 µl of the sample is extracted and the concentration taken from the standard curve is 0.15 ng/ml Noradrenaline.

Correction factor = 20/750 = 0.027

Concentration of the sample = 0.15 ng/ml x 0.027 = 0.004 ng/ml = 4 pg/ml Noradrenaline

### Conversion

Noradrenaline (ng/ml) x 5.91 = Noradrenaline (nmol/l)

Dopamine (ng/ml) x 6.53 = Dopamine (nmol/l)

## 7.1 Quality control

The confidence limits of the kit controls are indicated on the QC-Report.

## 8. Assay characteristics

Analytical Specificity (Cross Reactivity)	Substance	Cross Reactivity (%)	
		Noradrenaline	Dopamine
Analytical Specificity (Cross Reactivity)	Derivatized Adrenaline	0.14	0.03
	Derivatized Noradrenaline	100	0.87
	Derivatized Dopamine	0.2	100
	Metanephrine	< 0.003	< 0.007
	Normetanephrine	0.48	0.008
	3-Methoxytyramine	< 0.003	0.55
	3-Methoxy-4-hydroxyphenylglycol	0.01	< 0.007
	Tyramine	< 0.003	0.13
	Phenylalanine, Caffeinic acid, L-Dopa, Homovanillic acid, Tyrosine, 3-Methoxy-4-hydroxymandelic acid	< 0.003	< 0.007

Sensitivity (Limit of Detection)	Noradrenaline	Dopamine
	0.1 ng/ml x C*	0.25 ng/ml x C*

C\* = Correction factor (refer to 7.)

Analytical Sensitivity (750 µl undiluted sample)	Noradrenaline	Dopamine
	2.6 pg/ml	6.6 pg/ml

Functional Sensitivity (750 µl undiluted sample)	Noradrenaline	Dopamine
	4 pg/ml	10 pg/ml







Precision				
Intra-Assay Human EDTA-Plasma				
	Sample	Mean ± 3 SD (pg/ml)	SD (pg/ml)	CV (%)
Dopamine	high	1438.6 ± 465.6	155.2	10.8
	medium	565.9 ± 246.3	82.1	14.5
	low	56.4 ± 36.3	12.1	21.5
Noradrenaline	high	1377.4 ± 483.6	161.2	11.7
	medium	502.6 ± 126.9	42.3	8.4
	low	32.7 ± 15.3	5.1	15.6
Intra-Assay Cell Culture Medium (RPMI)				
	Sample	Mean ± 3 SD (pg/ml)	SD (pg/ml)	CV (%)
Dopamine	high	2784.5 ± 1238.7	412.9	14.8
	medium	1003.7 ± 526.2	175.4	17.5
	low	74.7 ± 51.6	17.2	23.0
Noradrenaline	high	2027.8 ± 712.5	237.5	11.7
	medium	716.5 ± 179.7	59.9	8.4
	low	46.0 ± 16.8	5.6	12.2

Recovery	Mean (%)	Range (%)	SD (%)	CV (%)
<b>Dopamine</b>				
Human EDTA-Plasma	97.7	83.7 – 115.9	11.8	12.1
Cell Culture Medium	98.6	77.7 – 113.4	12.1	12.2
<b>Noradrenaline</b>				
Human EDTA-Plasma	116.5	104.8 – 125.6	8.0	6.9
Cell Culture Medium	96.7	70.6 – 124.7	17.1	17.7

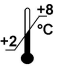

- **For literature or any other information please contact your local supplier.**
- **The liability of the manufacturer shall be limited to the replacement of defective products. The manufacturer takes no liability for any damages or expenses arising directly or indirectly from the use of this product.**

**Symbols:**

	Storage temperature		Manufacturer		Contains sufficient for <n> tests
	Expiry date	<b>LOT</b>	Batch code		
	Consult instructions for use	<b>CONT</b>	Content		
	Caution	<b>REF</b>	Catalogue number	<b>RUO</b>	For research use only!